

RescueLink® User's Guide

70-02084-31 A

Copyright© 2020 Cardiac Science Corporation. All rights reserved.



AT THE HEART OF SAVING
LIVES®

Information in this document is subject to change without notice. Names and data used in the examples are fictitious unless otherwise noted.

Trademark Information

Cardiac Science, the Shielded Heart logo, Powerheart, Mastertrak, MDLink, STAR, IntelliSense, Rescue Ready, RescueCoach, Rescuelink, and RHYTHMx are trademarks or registered trademarks of Cardiac Science Corporation. All other product and company names are trademarks or registered trademarks of their respective companies.

Copyright © 2020 Cardiac Science Corporation. All Rights Reserved.

Patents

This device is covered by the following U.S. and foreign patents:

5,792,190; 5,999,493; 5,402,884; 5,579,919; 5,749,902; 5,645,571;
6,029,085; 5,984,102; 5,919,212; 5,891,172; 5,674,266; 5,700,281;
5,891,173; 5,968,080; 6,263,239; 5,797,969; D402,758; D405,754;
5,909,138; 6,173,203; 6,088,616; 5,897,576; 5,955,956; 6,083,246;
6,064,909; 6,038,473; 5,868,794; 6,115,638; 6,366,809; 5,474,574;
6,246,907; 6,289,243; 6,411,846; 6,480,734; 6,658,290; EP00756878

Other U.S. and foreign patents pending.



Cardiac Science Corporation
500 Burdick Parkway
Deerfield, WI 53531 USA
800.426.0337 (USA)
262.953.3499
care@cardiacscience.com
www.cardiacscience.com

Contents

Getting Started.....	4
Communications Port Setup	5
Setting RescueLink Options.....	10
Setting Up Playback Controls	14
Using RescueLink	15
Clearing Rescue Data	15
Formatting Date and Time.....	16
Entering Annotation Data	16
Exporting RescueLink Files.....	18
Finding Rescue Data.....	19
Getting AED Date and Time.....	20
Getting Rescue Data.....	20
Opening Existing RescueLink Files.....	21
Playing Back Rescue Data.....	21
Printing Rescue Data	22
Setting AED Date and Time	24
Viewing Resue Data	25
Getting Help.....	31
Troubleshooting the AED.....	31
Troubleshooting RescueLink.....	33

Getting Started

The following drawings show the AED parts and their locations for the G3, G3 Plus, G3 Elite, and G3 Pro.

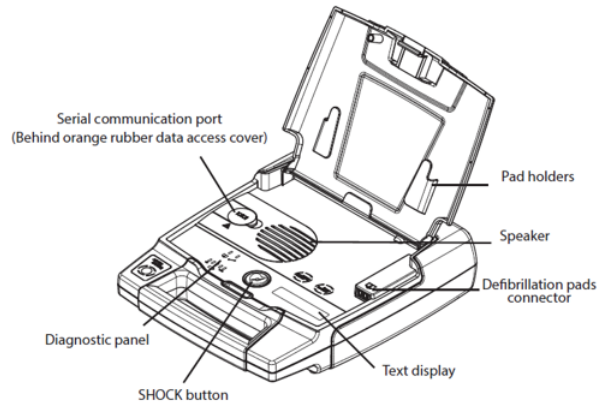


Figure 1: G3, G3 Plus, and G3 Elite

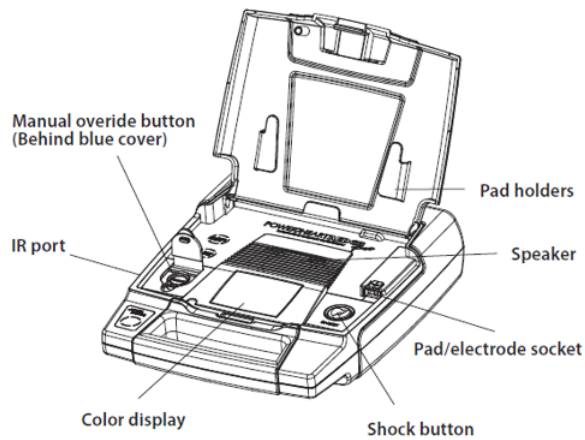


Figure 2: G3 Pro

Communications Port Setup

In order to download data from an AED, the communications port must be set up correctly. RescueLink® will automatically detect the communications port or it will allow the user to manually select the communications port.

Connecting the AED to your computer

In order for the rescue data to be transferred from the AED to a computer, the AED must be connected to your computer via the serial communication cable.

Note: If not already started, launch RescueLink.

For the G3 Pro:

- ◆ Connect the IrDa dongle into the IR port.


For the G3, G3 Plus, and G3 Elite:

- ◆ Connect one end of the serial communication cable to the serial communication port of the AED.
 - The serial connector of the Powerheart® G3 AED is located at the left side under the cover.
- ◆ Connect the other end of the serial communication cable to the appropriate communications port on the back of the computer.
 - The serial communication cable has a 9-pin connector. An adapter must be used if the communications port requires a 25-pin connector.


Automatically Detecting Communications Port

To have RescueLink automatically detect the communications port on the computer:

1. Connect the AED to your computer.
2. If not already started, launch RescueLink.
3. Open the AED Lid.
4. Select Auto Detect Communications Port from the Communications menu.

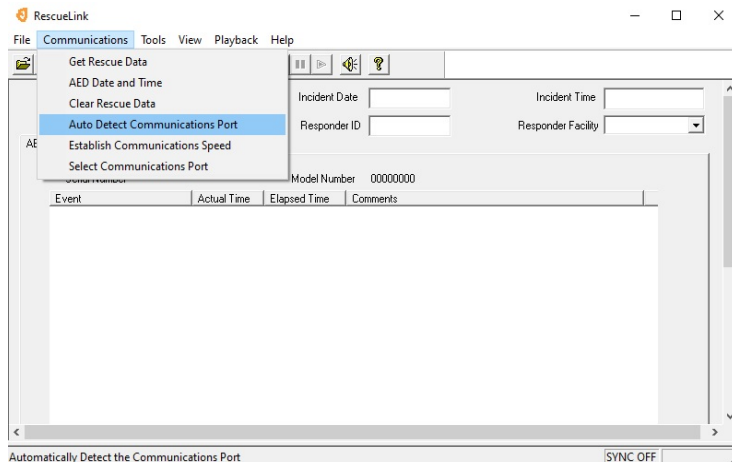


RescueLink will display the communications ports that are available, not available, and which communications port the AED is currently using.

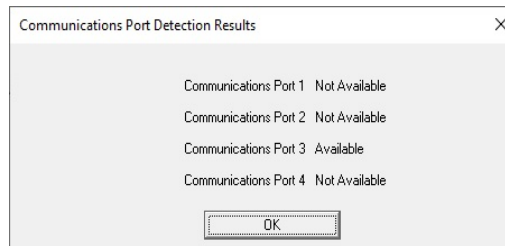


Manually Selecting Communications Port

To manually select a communications port on the computer:

1. Click on Select Communications Port from the Communications menu.

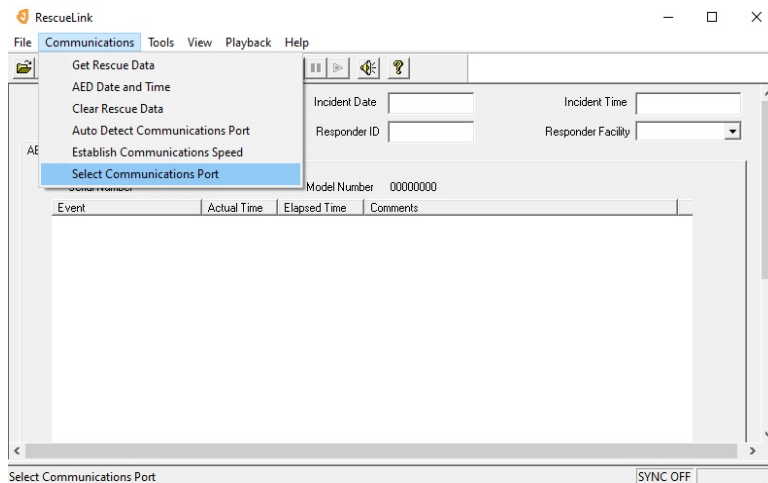


2. Select the desired communications port and click OK.

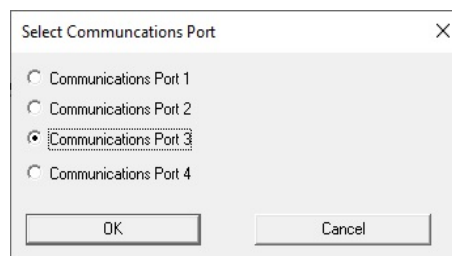


To make COM 5 to COM 8 available:

1. Click on the Tools menu, Options, check More Communications Ports, and click OK.

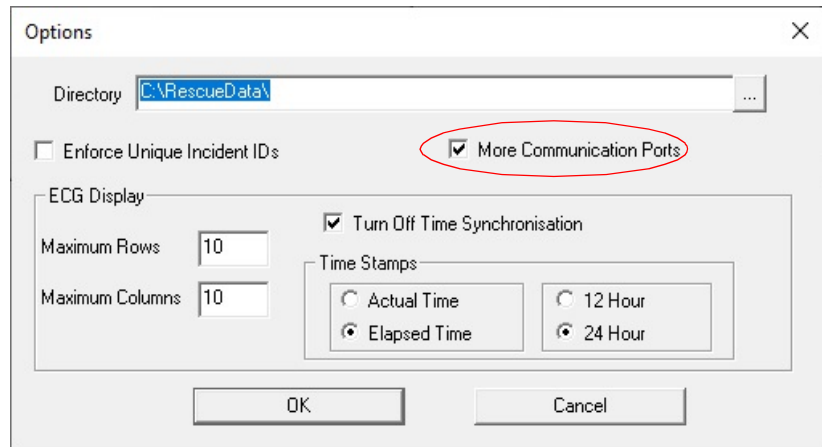

2. Click on Select Communication Port from the Communications menu. COM 1 to COM 8 display.
3. Click on the desired communications port and click OK.

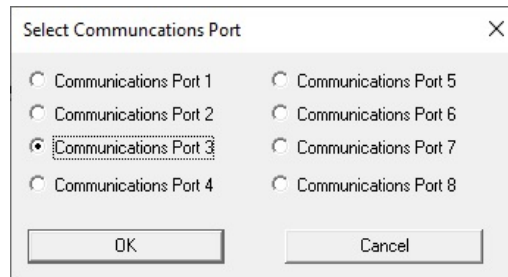


Note: Repeat steps 2 and 3 to select a different communications port.

Establishing Communications Speed

To establish communications between the AED and the computer:

1. Connect the AED to your computer.
2. Open the AED lid
3. Automatically or manually detect the communications port.
4. Click on Establish Communications Speed from the Communications menu.

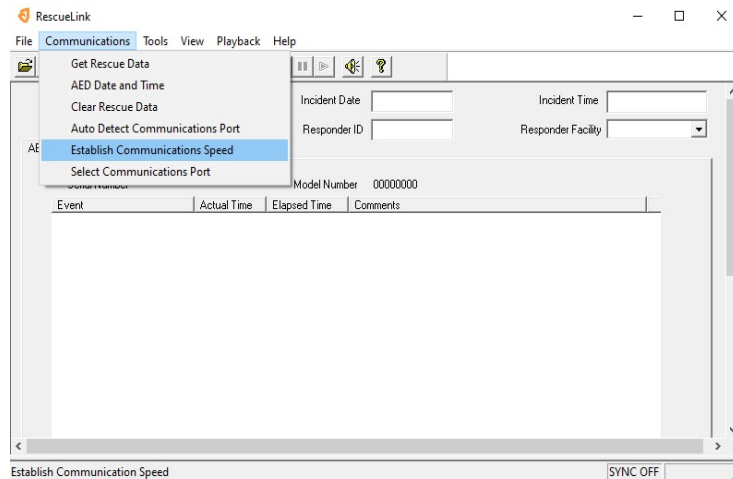

Setting RescueLink Options

ECG Display

RescueLink will allow the user to change the ECG rows and columns that are displayed on screen. This setting is also used for printing the ECG Heart Rhythm report.

1. Select Options from the Tools menu.
2. Enter a number 3 to 10 in the Maximum rows.

Note: The program defaults to 4.

3. Enter a number 3 to 30 in the Maximum columns.

Note: The program defaults to 10.

Data Directory

During installation, RescueLink installs a file called rescuelink.mdb in the following default location:


C:\RescueData\rescuelink.mdb.

The data directory or folder where the rescuelink.mdb file resides stores all ECG and incident data. In order for the user to view the ECG and incident data, the data directory field must point to the folder where the rescuelink.mdb resides. To change the data directory follow the instructions below.

Note: If the user selects a different path during installation, then the directory must point to the location of the rescuelink.mdb file.

To change the Data Directory location:

1. Verify ODBC data source has the desired location.
2. Select Options from the Tools menu.
3. The Directory displays the folder where the rescuelink.mdb resides. To change the location, select the icon at the end of the path to select the location.
4. Click OK.

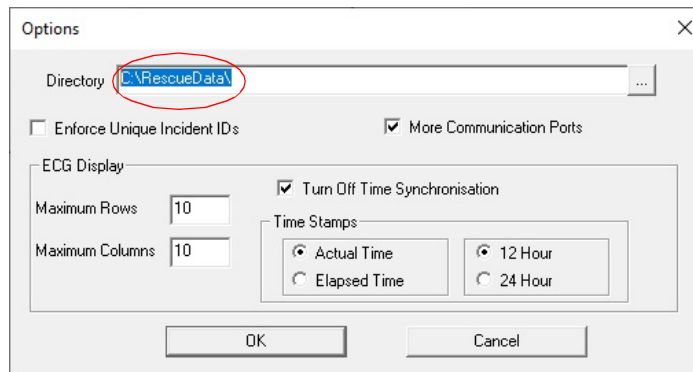

Note: RescueLink must be closed and then re-opened for the settings to take effect.


Turning On/Off Time Synchronization

RescueLink allows for synchronizing rescue data with the AED clock, provided the internal clock of the computer has been synchronized to the AED clock prior to the data transfer.

If time synchronization is enabled, then rescue data received from internal memory will be synchronized to the computer time at which the data transfer was initiated.

RescueLink will also display to the user if time synchronization is on or off in the status bar.

Note: The program installs with this option turned **off**.

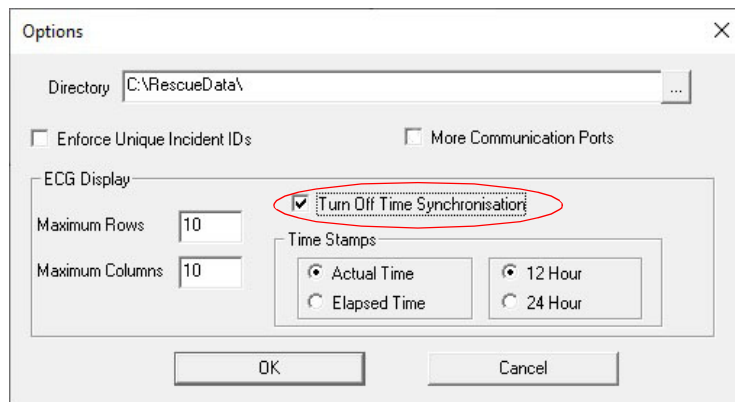

To turn Time Synchronization on:

1. Select Options from the Tools menu.
2. Un-check Turn off Time Synchronization.


To turn Time Synchronization off:

1. Select Options from the Tools menu.
2. Check Turn off Time Synchronization.

Formatting Time Stamps

To format time stamps:

1. Select Options from the Tools menu.
2. Select Actual Time or Elapsed Time.
3. Select 12 Hour or 24 Hour.

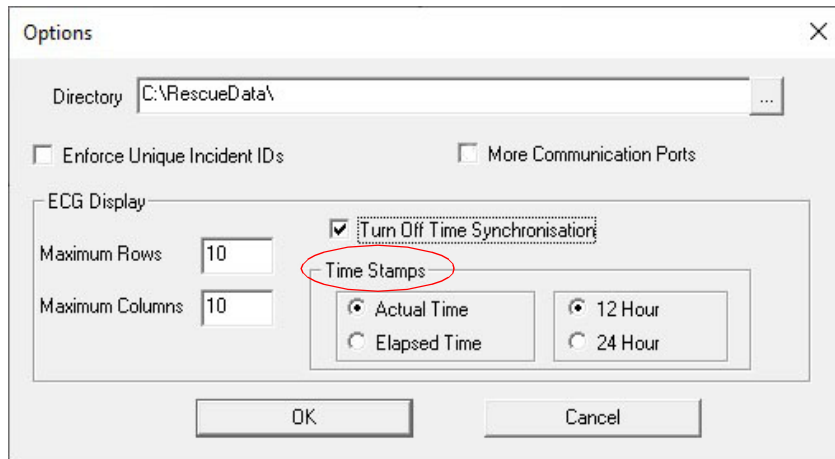


Setting Up Playback Controls

To set up or change the playback controls:

1. Select Playback Controls from the Playback menu.
2. Select the Scrolling Updates.

The scrolling update options are Fast, Medium, Slow, or None. The Scrolling update option is dependent upon the processing power of the computer.

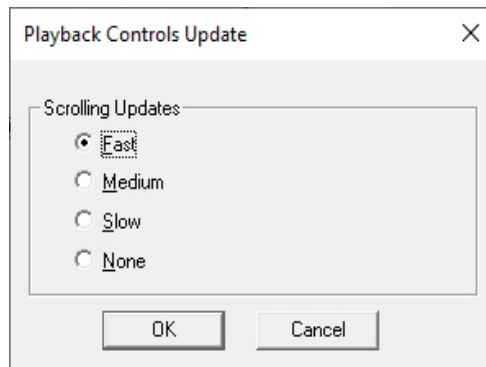

Using RescueLink

Clearing Rescue Data

Once the rescue data has been transferred from the AED to RescueLink, the data should be cleared from the AED's internal memory.

Clearing a Rescue From the Internal Memory of the AED

To clear rescue data from internal memory of the AED:

1. Connect the AED to your computer.
2. Open the AED lid.
3. Select Clear Rescue Data from the Communications menu.

Deleting Current Incident from RescueLink



CAUTION: A deleted incident is not recoverable.

- ◆ Select Delete Current Incident from RescueLink from the Tools menu

OR

- ◆ Select the Delete Rescue icon on the toolbar: 

Formatting Date and Time

RescueLink uses the short date format from Windows Control Panel Format.

RescueLink uses the time format from Regional Settings in the Windows Control Panel.

Note: A date or time field defaults to the current date and time of the computer.

Entering Annotation Data

RescueLink will allow the user to add, edit or delete user annotation data. It will not allow you to modify annotation recorded by the AED.

User annotation characteristics:

- ◆ User annotation will appear in the Event List Summary Report.
- ◆ User annotation will appear on the AED Data Tab in the Event List.
- ◆ The user can “drag and drop” the user annotation to any location on the ECG Heart Rhythm.
- ◆ The user can double click on the user annotation in the Event List and the event will appear in the upper left hand corner of the ECG Heart Rhythm.

Adding User Annotation

To add annotation data:

1. Right click on ECG Heart Rhythm area.
2. Select Add Annotation.
3. Enter Annotation data in text box.
4. Click OK.

The screenshot shows the RescueLink software interface. At the top, there is a menu bar with 'File', 'Communications', 'Tools', 'View', 'Playback', and 'Help'. Below the menu bar is a toolbar with various icons. The main window displays the 'AED Data' section with the following information:

Incident ID: sample1
Incident Date: 2/21/2007
Incident Time: 10:46:24 AM
Responder ID:
Responder Facility:

AED Data

Serial Number	Model Number
00361542	00019300

Event	Actual Time	Elapsed Time	Comments
Lid Open	10:46:24 AM	00:00:00	
Electrodes Placed	10:46:29 AM	00:00:05	
Start of Analysis	10:46:34 AM	00:00:10	Analysis #0
Start of Charge	10:46:39 AM	00:00:15	
Shock Advised	10:46:39 AM	00:00:15	
End of Charge	10:46:45 AM	00:00:21	
Shock Synchronised	10:46:50 AM	00:00:26	1 (Biphasic) 194 joules 53 ohms
Start CPR	10:47:00 AM	00:00:36	
End CPR	10:49:00 AM	00:02:36	
Start of Analysis	10:49:07 AM	00:02:43	Analysis #1
Start of Charge	10:49:11 AM	00:02:47	
Shock Advised	10:49:11 AM	00:02:47	
End of Charge	10:49:19 AM	00:02:55	
Shock Synchronised	10:49:24 AM	00:03:00	2 (Biphasic) 259 joules 53 ohms

Below the table is a grid of ECG waveforms. A context menu is open over the grid, with the following options: Play, Playback Controls, Invert Rhythm, Time Stamps, Event Markers, Add Annotation (highlighted), Edit Annotation, and Properties.

The screenshot shows the 'Annotation Data' dialog box. It has a title bar with a close button (X). Inside the dialog, there is a text input field. Below the input field are three buttons: 'Delete', 'OK', and 'Cancel'.

Editing User Annotation

To edit user annotation data:

1. Right click on user annotation data.
2. Select Edit Annotation.
3. Edit Annotation data in text box.
4. Click OK.


Deleting User Annotation

To delete user annotation data:

1. Right click on user annotation data.
2. Select Edit Annotation.
3. Click Delete button.
4. Click OK.

Exporting RescueLink Files

To export or save a file to a .svl file format:


- ◆ Select Export (.svl) file from File menu.
OR select the Export Rescue option from the toolbar: 
- ◆ Enter the desired filename and location to save the file.
- ◆ Click Save.

Finding Rescue Data

To find specific data:

1. Select Find Incident from Tools menu.
2. Select the one of the following options:
 - Incident ID
 - Incident Date
3. Click (or double-click) on an item to highlight the desired data and click OK.

RescueLink will display the specified incident.

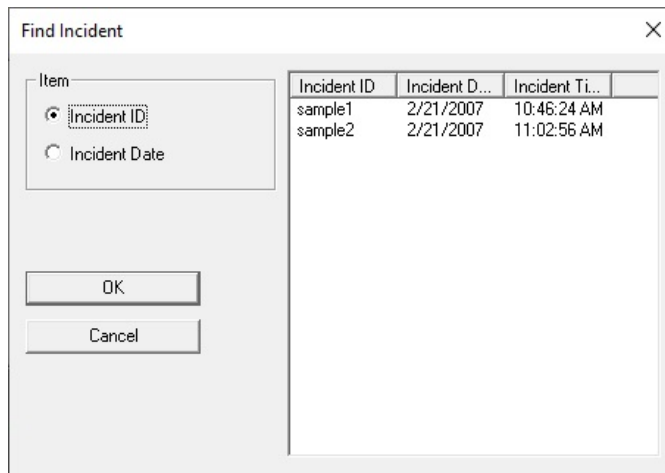


Getting AED Date and Time

1. Connect the AED to your computer.
2. Open the AED Lid.
3. Select AED Date and Time from the Communications menu.
4. Click on Get button.

RescueLink will display the date and time that is currently stored in the AED.

Getting Rescue Data

Rescue data can be transferred from a Cardiac Science® AED and stored in an incident file.


1. Connect the AED to your computer.
2. Open the AED lid.
3. Select Get Rescue Data from the Communications menu.
OR right click on the ECG graph area and select Get Rescue Data.
OR click the Get Rescue Data tool bar icon: 

Opening Existing RescueLink Files

The user can open existing RescueLink files created in earlier versions of RescueLink.

RescueLink will open a file with (.svl) extension and merge it into the new format for the latest version of RescueLink to read. The old file will remain on the computer.


To open an existing RescueLink file:

1. Select Open (.svl) file from File menu.
OR select the Open (.svl) file option from the toolbar: 
2. Select the desired filename and location of the (.svl) file.
3. Click Open.


Playing Back Rescue Data

RescueLink will play, pause, or stop ECG data of a rescue.


To play:

- ◆ Select Play from the Playback menu.
OR use the mouse to right click on the ECG area and select Play.
OR use the Play icon on the toolbar: 

To pause:

- ◆ Select Pause from the Playback menu.
OR use the mouse to right click on the ECG area and select Pause.
OR use the Pause icon on the toolbar: 

To stop:

- ◆ Select Stop from the Playback menu.
OR use the mouse to right click on the ECG area and select Stop.
OR use the Stop icon on the toolbar: 


Printing Rescue Data

RescueLink allows the user to print four predefined reports.

ECG Heart Rhythm

Note: If Time Stamps and Event Markers are not currently displayed on the screen, they will not be printed on the report.

To print the ECG Heart Rhythm:

1. Select Print Reports from the File menu.
OR Select the Print Reports toolbar icon: 
2. Click on the Report Setup.
3. Select ECG Heart Rhythm.
4. Click OK.
5. Click Print.

Print Setup

RescueLink will allow the user to set up the print options for the ECG Heart Rhythm report **only**. All other reports generated by RescueLink have a predefined paper size and orientation setting.


To setup print options for the ECG Heart Rhythm only:

- ◆ Select Print Setup from the File menu.

Event List Summary

Note: The report will print in Portrait Mode only.


To print the Event List Summary:

1. Select Print Reports from the File menu.
OR Select the Print Reports toolbar icon: 
2. Click on the Report Setup.
3. Select Event List Summary.
4. Click OK.
5. Click Print.

AED Operating Parameters

Note: The report will print in Portrait Mode only.


To print AED Operating Parameters:

1. Select Print Reports from the File menu.
OR Select the Print Reports toolbar icon: 
2. Click on the Report Setup.
3. Select AED Operating Parameters.
4. Click OK.
5. Click Print.

Battery Information

Note: The report will print in Portrait Mode only.


To print Battery Information:

1. Select Print Reports from the File menu.
OR Select the Print Reports toolbar icon: 
2. Click on the Report Setup.
3. Select Battery Information.
4. Click OK.
5. Click Print.

Setting AED Date and Time

To set AED Date and Time:

1. Connect the AED to your computer.
2. Open the AED lid.
3. Select AED Date and Time from the Communications menu.
4. Click on Set button.
5. The AED will be set to the computer's current date and time.

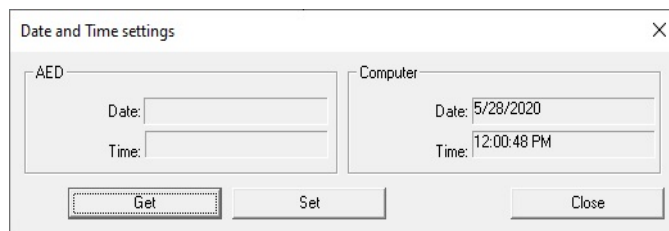

Viewing Resue Data

RescueLink allows the user to view rescue data downloaded from a Cardiac Science AED.

ECG Heart Rhythm

RescueLink displays the ECG heart rhythm recorded by Cardiac Science AEDs. The ECG Heart Rhythm is displayed on the main screen below the splitter. The user can drag the horizontal splitter to fully display or hide the ECG heart rhythm.

Important: Rescue data managed by RescueLink is for archival purposes only. RescueLink does not attempt to interpret medical information and is not a medical device.

The ECG Heart Rhythm contains the following properties:

- ◆ The ECG Grid displays 1 second increments.
- ◆ Events are marked with arrows and a text description.
- ◆ A scroll bar and arrow key allow the user to browse the rescue data.
- ◆ Double clicking on an event in the AED Data tab will display the specific event in the upper right corner of the ECG grid area.

AED Operating Parameters

The user will be able to print or preview the following AED Operating Parameters for the current incident displayed on the screen:

- ◆ Model Number
- ◆ Serial Number
- ◆ Code Version
- ◆ Energy Protocol
- ◆ Maximum shocks per rescue
- ◆ Same Energy After Conversion
- ◆ Daylight Saving Time
- ◆ Electrode Test
- ◆ Prompt Sound Volume
- ◆ SVT Rate
- ◆ VF/VT Rate
- ◆ Continue CPR Prompt
- ◆ CPR Time
- ◆ User Defined Identification
- ◆ Initial Prompt (G3 Plus and G3 Elite)
- ◆ Compression-only CPR (G3 Plus and G3 Elite)
- ◆ CPR First (G3 Plus and G3 Elite)
- ◆ Manual Button (G3 Pro)
- ◆ Remain in Manual Button (G3 Pro)

Preview AED Operating Parameters

To preview AED Operating Parameters:

- ◆ Select AED Operating Parameters from the Tools menu.

Battery Information

The user can view or print the following battery information for the current incident displayed or directly from the AED:

- ◆ Battery Type
- ◆ Date of Installation
- ◆ Number of Charges Completed
- ◆ Time of Operation
- ◆ Days of Standby Operation
- ◆ Battery Capacity Remaining
- ◆ Battery Identification if available

Preview Battery Information

To preview Battery Information:

1. Select Battery Information from the Tools menu.
2. Select Current Incident or AED.

Events List Summary

The user will be able to print or preview the event list for the current incident displayed on the screen.

Preview the Event List Summary

To preview the Event List Summary:

- ◆ Click on the AED Data tab.

Self-test History Information

The user will be able to print or preview the AED Self-test History for the current incident displayed on the screen or from an AED. It will display the last 20 self-test events from the AED.

Preview AED Self-Test History

To preview the AED Self-test History:

1. Select AED Self-test History from the Tools menu.
2. Select Current Incident or AED.

Event Markers

RescueLink allows the user to turn on or off the event markers.

To turn on or off the event markers:

- ◆ Select Event Markers from the View menu.
OR right click on the ECG area and select Event Markers.

Note: If Event Markers is checked, they will be displayed on the screen and also on the ECG Heart Rhythm report.

Time Stamps

The RescueLink displays time stamps in one second intervals from the AED lid open time to lid closed.

RescueLink allows the user to turn on or off the time stamps.

To turn on or off the time stamps:

- ◆ Select Time Stamps from the View menu.
OR right click on the ECG area and select Time Stamps.

Note: If Time Stamps is checked, they will be displayed on the screen and on the ECG Heart Rhythm report.

Invert Rhythm

The heart rhythm displayed in RescueLink may appear inverted. RescueLink will allow the user to view the heart rhythm inverted.

To invert the ECG heart rhythm:

- ◆ Select Invert Rhythm from the View menu.

OR right click on the ECG area and select Invert Rhythm.

Note: If Invert Rhythm is checked, the heart rhythm will be displayed inverted and also appear inverted in the ECG Heart Rhythm report.

Status Bar

RescueLink allows the user to turn on or off the status bar.

To turn on or off the status bar:

- ◆ Select Status Bar from the View menu.

Note: If Status Bar is checked, it will be displayed on the screen.

Toolbar

RescueLink allows the user to turn on or off the toolbar. The toolbar contains buttons that provide short cuts to some of the common commands.

To turn on or off the toolbar:

- ◆ Select Toolbar from the View menu.

Note: If Toolbar is checked, it will be displayed on the screen.

ECG Properties

RescueLink allows the user to view the ECG Properties.

This includes the following:

- ◆ Length of recording (in seconds)
- ◆ Sound recorded
- ◆ Number of annotations
- ◆ Number of shocks

To view ECG properties:

- ◆ Right click on the ECG area
- ◆ Select ECG Properties

Getting Help

Troubleshooting the AED

AED Error Codes

During its normal use or self-testing, the AED may detect errors and record them in the AED Self-Test History. The AED error codes and their expected actions are listed below. If an error code appears on your self-test history that is not on the list, contact Technical Support.

Battery Alert (0x32, 0x33, 0x39)

The battery in the AED needs to be replaced.

Refer to the user manual of the AED for instructions on replacing the battery.

AED Operating Temperature Alert (0x40, 0xE1)

Exposing the Powerheart AED with the battery installed to temperature extremes outside the normal operating and standby conditions will cause the self-tests to be disabled and could cause the AED to function improperly. Storing the AED outside these conditions for 5 consecutive days will result in a "Service" alert.

Operating Temperature of AED Without Display -40°C to 65°C (-44°F to 149°F)

Operating Temperature of AED With Display -30°C to 65°C (-22°F to 149°F)

Electrodes Alert (0x55)

The electrodes are either missing from the AED or need to be replaced.

Refer to the user manual of the AED for instructions on installing electrodes in the AED.

AED Alert (0xD0)

The AED has detected an error. Contact Technical Support.

AED Display Alert (0x5A)

The AED has detected an error with the display. Contact Technical Support.

Shock Button Alert (0x71)

The rescue button was depressed while the AED was performing self-testing.

Close and open the lid of the AED, ensuring the rescue button is not depressed until the LEDs on the diagnostic panel stop blinking.

AED Lid Alert (0xE0)

The lid of the AED was not completely closed. Ensure the AED lid can be closed completely when the AED is not in use.

Serial Communication Cable Alert (0x5F)

The serial communication cable was connected to the serial communication port during self-testing of the AED. The cable should be disconnected from the AED when not being used to retrieve rescue data.

Troubleshooting RescueLink

If you receive an Error Code while using RescueLink, see the corresponding section below.

Communications Error (1, 2, 4, 13, 14, 15, 16, 17, 18, 19, 21, 26, 54, 66, 67, 68, 94)

The following reasons could cause communication problems between the AED and the computer:

- ◆ Computer's communications port is configured incorrectly.
 1. Use RescueLink to correctly set up communications port.
 2. Make sure the computer's communications port is correctly configured.
- ◆ AED is not in Communications Mode.
- ◆ Wrong serial cable is connected.

If these solutions do not work, contact Technical Support for assistance.

Print Error (32, 56, 57, 65)

The following reasons could cause printing problems:

- ◆ ODBC and data directory are not pointing to the same location.
- ◆ Report is missing from directory.
- ◆ Printer is not installed on computer or is configured incorrectly.

Import Error (8, 24, 52, 134, 135, 136, 137, 139, 140, 141, 142, 143, 144, 145, 146, 147, 148, 149, 150, 151, 152, 153, 154, 155, 156, 157)

The following reasons could cause importing problems:

- ◆ ODBC and data directory are not pointing to the same location.

Write Error (3, 23, 27, 30, 35, 36, 37, 38, 39, 42, 44, 47, 48, 51, 74, 75, 79, 80, 90, 91, 92, 97, 98, 99, 100, 101, 102, 103, 104, 105, 106, 107, 108, 109, 110, 111, 112, 113, 114, 115, 116, 117, 118, 119, 120, 121, 122, 123, 124, 125, 126, 127, 128, 129, 130, 131, 132, 133, 138)

You do not have write access to the RescueLink folder. Contact your network administrator to obtain write access privileges.

If this solution does not work, contact Technical Support for assistance.

ODBC Error (33)

ODBC and Data Directory are not pointing to same rescuelink.mdb file.

Verify ODBC is set up correctly.

Verify Data Directory is pointing to correct file.

If this solution does not work, contact Technical Support for assistance.

Missing Report Error (56, 57)

- ◆ Verify the following reports exist on the computer. If they do not exist, reinstall RescueLink:
 - svlrpt1.dsr
 - svlrpt2.dsr
 - svlrpt3.dsr
- ◆ Verify the aforementioned files reside in the same folder as the rescuelink.mdb file.

If this solution does not work, contact Technical Support for assistance.

Database Error (63, 159)

- ◆ The rescuelink.mdb file is not located in the directory as specified in ODBC and the data directory. Verify that the file exists. If the file does not exist, contact Technical Support. It is recommended that you backup a copy of this file on a regular basis.
- ◆ Verify ODBC is set up correctly.
- ◆ Verify Data directory is pointing to correct file.

If these solutions do not work, contact Technical Support for assistance.

Communications Mode Error (28)

AED is not in Communications Mode

- ◆ Connect the AED to your computer.
- ◆ Open the AED lid.

If this solution does not work, contact Technical Support for assistance.

Missing ECG File Error (5, 9 , 10, 11, 12, 45, 58, 95, 161)

- ◆ The incident is not in the data directory.
- ◆ The incident is in use by another user.
- ◆ ODBC and Data Directory are not pointing to same rescuelink.mdb file.

If this solution does not work, contact Technical Support for assistance.

Export Error (53, 96)

The incident is missing the ECG file or it is in use by another user.

Read Only Error (60, 69, 162)

The file rescuelink.mdb is marked as Read-only. Locate the file using Windows Explorer and uncheck the Read-only property of this file.

If this solution does not work, contact Technical Support for assistance.

Memory Error (76, 86)

Unable to allocate enough memory to finish the operation.

Try closing other programs that are currently running on your computer.

If this solution does not work, contact Technical Support for assistance.

Option Unavailable (55, 73)

The AED is not equipped with this option.

Not Used (20, 22, 25, 59, 61, 62, 71, 160)

This error code is not currently used in this program.

Technical Support (0, 29, 40, 41, 64, 65, 70)

This issue requires calling Technical Support to resolve.

User Annotation Error (93)

Use your mouse to right click within 1 second of the user annotation marker.

Edit Error (43, 49, 50, 158, 163, 164, 165)

If you are unable to edit while performing an operation, then wait until the operation is completed.

Cardiac Science Corporation • 500 Burdick Parkway, Deerfield, WI 53531 USA • 262.953.3500 • US toll-free 800.426.0337 •
Fax: 262.953.3499 • care@cardiacscience.com

Orders and Customer Care (US and International) • US toll-free 800.426.0337 • Fax: 262.953.3499 • care@cardiacscience.com

Technical Support • US toll-free 800.426.0337 • (US) Fax: 262.798.5236 • techsupport@cardiacscience.com •
(International) internationalservice@cardiacscience.com

Cardiac Science, the Shielded Heart logo, Powerheart, Mastertrak, MDLink, STAR, IntelliSense, Rescue Ready, RescueCoach, Rescuelink, RHYTHMx, are trademarks of Cardiac Science Corporation. Copyright © 2020 Cardiac Science Corporation. All Rights Reserved.



70-02084-31 A

